



Product catalogue

SmartIntegrate

Integrative Organisational Solutions for Laboratory Workstations



Table of Contents

Your Advantage in Everyday Laboratory Work.....	3
SmartIntegrate Integrative Organisational Solutions.....	6
The Configuration.....	6
Waldner SmartRail (work benches).....	7
Waldner SmartEdge (fume hoods).....	10
Waldner SmartFrame (fume hoods).....	28
Universal solutions (work benches/fume hoods).....	30
Additional Accessories.....	34
Customising.....	35
Item Overview.....	37
Technical Information.....	40

Your Advantage in Everyday Laboratory Work

Better Basics Laborbedarf – Your Advantage in Everyday Laboratory Work

Better Basics Laborbedarf optimises laboratories and their workspaces with innovative, award-winning organisational systems. We support our customers in creating the best possible working environment for their laboratory staff, helping to unlock the full potential of their laboratories. This approach ensures consistently high quality, efficiency, and process stability in the lab while providing a safe, ergonomic workplace with satisfied employees.

Enhancing Employee Satisfaction

Our SmartRack® and SmartIntegrate product families offer a comprehensive system to optimise laboratory workflows. With over 100 modular solutions, they not only ensure order but also noticeably improve working conditions. Enhanced ergonomics, reduced movement paths, and clearly structured workstations increase both efficiency and employee satisfaction. Safe and organised environments reduce stress and promote job satisfaction. Furthermore, they demonstrably contribute to improved workplace safety and occupational health through thoughtful placement and ease of use.


A System for Sustainable Process Improvement

Our product families are designed to standardise workflows and save time. The modular concept of SmartRack® and SmartIntegrate allows for flexible adaptations to any laboratory environment – from safety workbenches to large open-plan labs. These systems not only optimise space utilisation but also increase efficiency by an average of 30% and save up to 80% space. With a wide range of modules and accessories, they offer customised solutions for every need.

Efficient, Safe, and Satisfying

For decision-makers and managers, SmartRack® and SmartIntegrate ensure process stability and economical resource use. At the same time, these systems provide innovation managers and organisational developers with the ideal basis for forward-looking workplace designs. For laboratory staff, this means fewer errors thanks to well-structured work areas, a more relaxed working environment, and often the ability to finish earlier due to significantly improved workflows.

Discover how Better Basics Laborbedarf can transform your everyday laboratory work – for more efficient processes, satisfied employees, and excellent results.

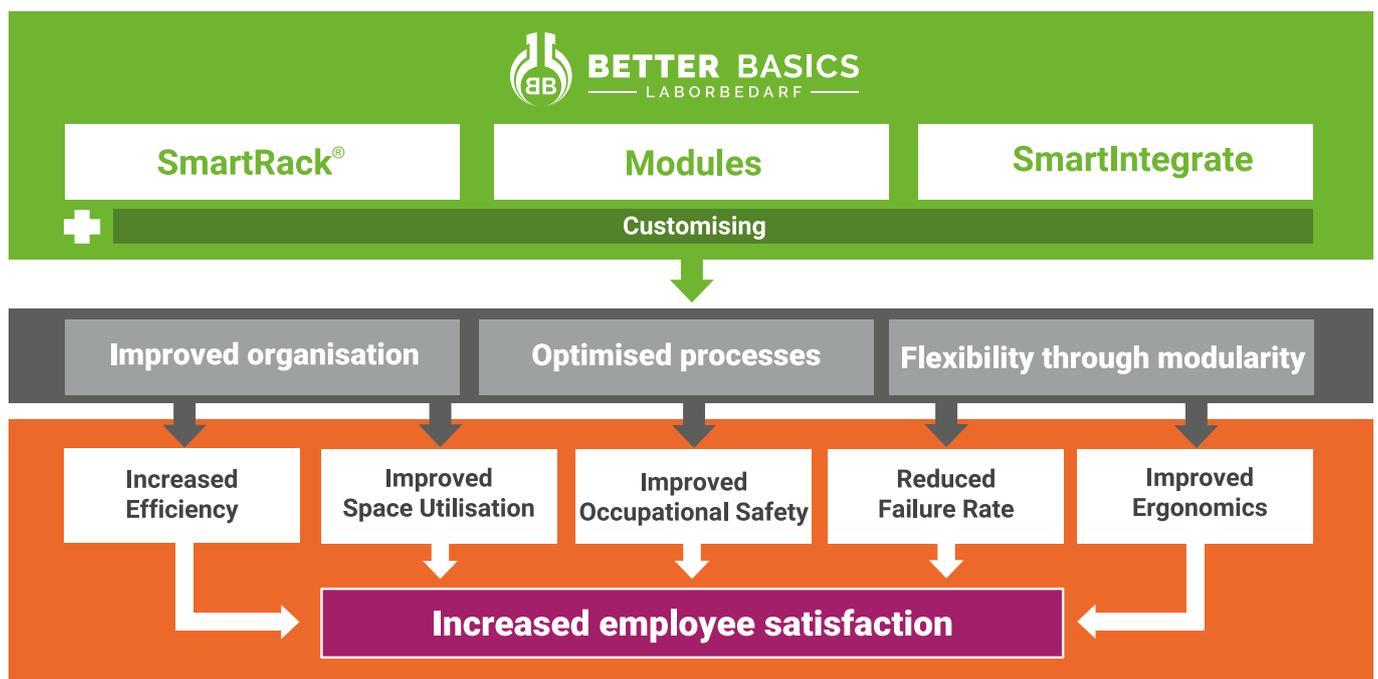

The product overview of Better Basics

▶ SmartRack®, the mobile solution for: laboratory workbench, fume hood and safety cabinet



+ Customised solutions

▶ SmartIntegrate, the integrative solutions for

I. Laboratory workbench



II. Fume hood



III. Safety cabinet



IV. Integration into laboratory equipment



+ Customised solutions

▶ Modules

- Currently, more than 100 modules
- All modules in the Better Basics catalogue can be used universally in all SmartRack® and SmartIntegrate solutions
- Customised modules can be developed based on the specific customer needs, all customised modules can be used in all SmartRack® and SmartIntegrate solutions



+ Customised solutions

SmartIntegrate - For Optimal Integration

The SmartIntegrate product line offers organisational solutions that integrate seamlessly into your existing lab infrastructure. The SmartFrame, for example, can be easily integrated into any fume hood, while the SmartRail provides additional storage space on media cells and laboratory benches. The SmartEdge optimises organisation in biosafety cabinets or fume hoods.

All SmartIntegrate organisational systems benefit from the compatibility of all Better Basics modules with these SmartIntegrate organisational solutions, ensuring intuitive and safe usage.

▶ **SmartIntegrate**, the integrative solutions for

I. Laboratory workbench



II. Fume hood



III. Safety cabinet



IV. Integration into laboratory equipment



✚ **Customised solutions**



Photo: The Waldner SmartRail WA 003, mounted in the integral rail of a Waldner SCALA media cell

SmartIntegrate | The Configuration

<p>1. Selection of workplace: Media cell/workbench, Fume hood</p>		
<p>2. Selection of the SmartIntegrate Version: SmartRail, SmartFrame, SmartEdge</p>		 
<p>3. Selection of the modules: Holders and storage box and storage box inserts (e.g. refer to catalogue „Modules“)</p>		 



Photo: SmartEdge WA 001 in a Waldner SCALA 1500 RWI fume hood



SmartRail WA 001




SmartRail WA 001 mounted, Demonstrators and modules are not included in the scope of delivery.

The SmartRail WA 001 has been specifically developed for use on WALDNER laboratory furniture, particularly for workbenches with shelves that have the WALDNER functional rail, such as for instance the WALDNER Simply Flex Bench. The product consists of a black anodised aluminium slot rail and side brackets made of PA12 and stainless steel, providing secure magnetic fixation. The SmartRail WA 001 can be additionally equipped with a SmartLock slide locking mechanism to prevent unintentional removal and ensure stable mounting on the WALDNER functional rail. The SmartLock is included in the scope of delivery. The versatile SmartRail offers flexible placement options in the laboratory. A larger number of SmartRails provides various module placement options and allows the laboratory environment to be tailored to specific requirements. All modules are fully compatible with this SmartRail.



SmartRail WA 001 mounted, without equipment.

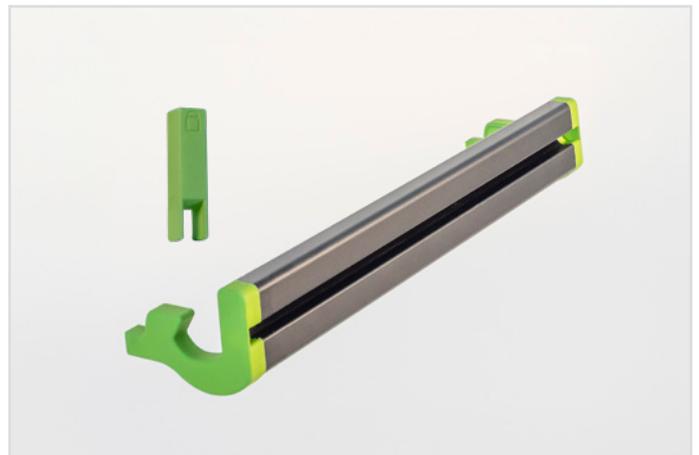

SmartRail WA 001, corner view.

Item	Type	PU	Item no.	Price in EUR
SmartRail WA 001	300 mm	1	03002	69,99
SmartRail WA 001	600 mm	1	03003	89,99
SmartRail WA 001	900 mm	1	03004	109,99
SmartRail WA 001	1200 mm	1	03005	129,99
SmartRail WA 001	1500 mm	1	03006	149,99



SmartRail WA 002



The SmartRail WA 002 has been specifically developed for use on WALDNER laboratory furniture, particularly for workbenches with shelves that have the WALDNER functional rail, such as for instance the WALDNER Simply Flex Bench. The product consists of a black anodised aluminium slot rail and side brackets made of PA12 and stainless steel, equipped with a robust twist lock that ensures stable mounting on the WALDNER functional rail. The versatile SmartRail offers flexible placement options in the laboratory. A larger number of SmartRail provides various module placement options and allows the laboratory environment to be tailored to specific requirements. All modules are fully compatible with this SmartRail.

SmartRail WA 002 mounted, Demonstrators and modules are not included in the scope of delivery.



SmartRail WA 002 mounted, without equipment.



SmartRail WA 002, corner view.

Item	Type	PU	Item no.	Price in EUR
SmartRail WA 002	300 mm	1	03007	79,99
SmartRail WA 002	600 mm	1	03008	99,99
SmartRail WA 002	900 mm	1	03009	119,99
SmartRail WA 002	1200 mm	1	03010	139,99
SmartRail WA 002	1500 mm	1	03011	159,99



SmartRail WA 003



SmartRail WA 003 mounted, Demonstrators and modules are not included in the scope of delivery.

The SmartRail WA 003 has been specifically developed for use with WALDNER laboratory furniture, particularly for workbenches with classical functional rail and shelves with classical functional rail of the SCALA service spine series. The product consists of a black anodised aluminium slot rail and side brackets made of PA12 and stainless steel, equipped with the original WALDNER clamp adapter to ensure stable mounting on the slot rail and shelves with the classical WALDNER functional rail. The versatile SmartRail offers flexible placement options in the laboratory. A larger number of SmartRail provides various module placement options and allows the laboratory environment to be tailored to specific requirements. All modules are fully compatible with this SmartRail.



SmartRail WA 003 mounted, without equipment.



SmartRail WA 003, corner view.

Item	Type	PU	Item no.	Price in EUR
SmartRail WA 003	300 mm	1	03012	79,99
SmartRail WA 003	600 mm	1	03013	99,99
SmartRail WA 003	900 mm	1	03014	119,99
SmartRail WA 003	1200 mm	1	03015	139,99
SmartRail WA 003	1500 mm	1	03016	159,99



SmartEdge elite WA 001



SmartEdge elite WA 001 in a WALDNER SCALA fume hood, Slot rails are not included in the scope of delivery.

The SmartEdge elite WA 001 is compatible with all of the following Waldner SCALA fume hoods:

- Bench-mounted fume hood with rear wall installation
- Secuflow bench-mounted fume hood with rear wall installation
- Low ceiling bench-mounted fume hood with rear wall installation
- Secuflow low ceiling bench-mounted fume hood with rear wall installation
- Secuflow bench-mounted fume hood with side installation
- Secuflow low ceiling bench-mounted fume hood with side installation
- Walk-in fume hood with side installation



SmartEdge elite WA 001 in a WALDNER SCALA fume hood, Slot rails, demonstrators and modules are not included in the scope of delivery.

The SmartEdge elite WA 001 is a modular organization system for installation in WALDNER laboratory fume hood of the SCALA series. It consists of two stainless steel L-profiles that are permanently connected to the stand mounts on the rear wall of the fume hood. The system allows flexible adaptation to the width and type of the respective fume hood.

The L-profiles can be mounted either on the outer stand rod holders on the left and right or, if available, on the middle rows of rod holders. This allows configuration across the entire width or only across parts of the fume hood rear wall.

Cross connectors in the form of precision-manufactured slot rails create a highly variable organizational surface on the rear wall of the fume hood, which can be adapted precisely to individual work requirements and the existing fume hood architecture. All modules are fully compatible with the slot rails for SmartEdge elite WA 001. (Slot rails and modules are not included in the scope of delivery.)

The SmartEdge elite WA 001 ensures stable mounting, high load capacity, and efficient use of space in WALDNER SCALA series laboratory fume hoods.

The functionality of the fume hood is not restricted by the installation of the SmartEdge elite WA 001.

No structural changes to the fume hood are required for the installation of the Smart Edge elite WA 001. After dismantling the SmartEdge elite WA 001, the fume hood is returned to its original condition.


Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99



Configuration and selection of your SmartEdge elite WA 001 for Waldner SCALA

1. Check whether your fume hood is equipped with a rear wall installation or a side installation.
2. Each laboratory fume hood has a specific width. You can find this width on the type plate of the fume hood. Typically (except in the case of customised designs), fume hoods have dimensions of 1200 mm, 1500 mm, 1800 mm, 2100 or 2400 mm.
3. The two stainless steel L-profiles can be used for **all** fume hood widths. Depending on the width of the fume hood and the desired position of the L-profiles, you can select the required number of matching slot rails.
4. Please note that the stainless steel L-profiles can be fixed in a total of 49 different mounting positions on the stand rod holders, depending on the fume hood width and the selected installation position.

Waldner SCALA fume hood with rear wall installation - Width 1200 mm

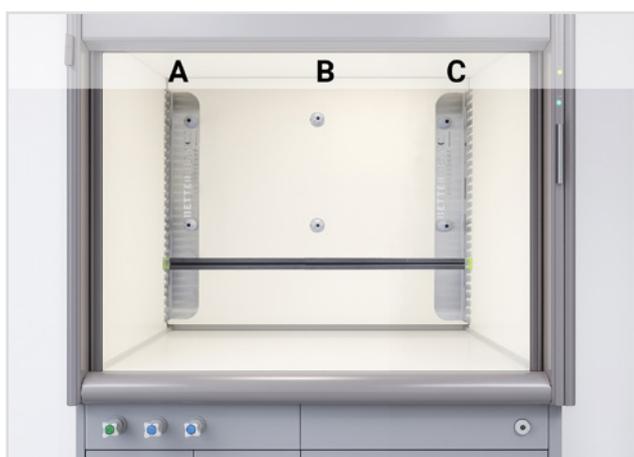


Fume hoods with a width of 1200 mm have three rows (A, B and C) of stand rod holders.

To connect the outer rows of stand rod holders (A-C), select the required number of slot rails with a width of 1101 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	1101 mm (A-C)	1	04320	69,99

Waldner SCALA fume hood with rear wall installation - Width 1200 mm

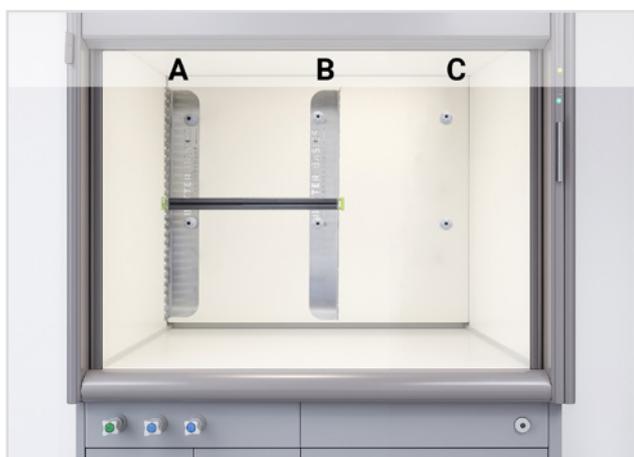

Fume hoods with a width of 1200 mm have three rows (A, B and C) of stand rod holders.


To connect the middle row of stand rod holders with one of the outer rows (A-B, B-C), select the required number of slot rails with a width of 613 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	613 mm (A-B, B-C)	1	04307	54,99



Waldner SCALA fume hood with side installation - Width 1200 mm

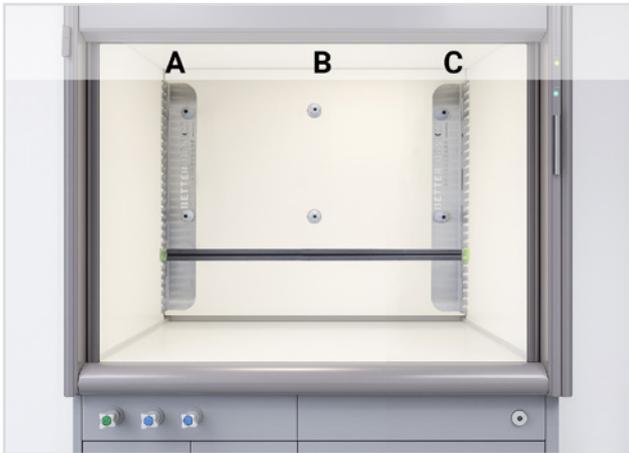


Fume hoods with a width of 1200 mm have three rows (A, B and C) of stand rod holders.

To connect the outer rows of stand rod holders (A-C), select the required number of slot rails with a width of 901 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	901 mm (A-C)	1	04318	64,99

Waldner SCALA fume hood with side installation - Width 1200 mm

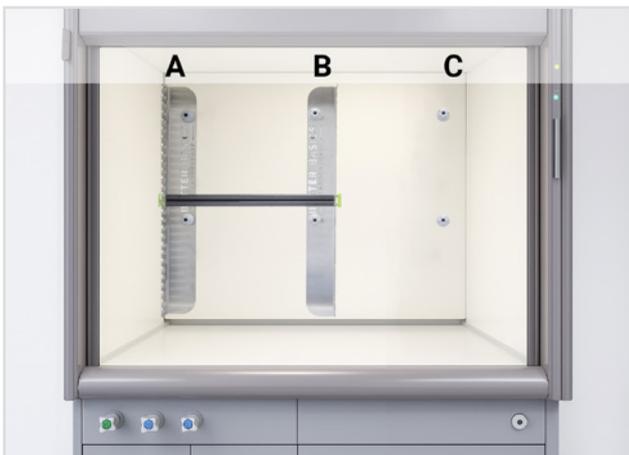

Fume hoods with a width of 1200 mm have three rows (A, B and C) of stand rod holders.


To connect the middle row of stand rod holders with one of the outer rows (A-B, B-C), select the required number of slot rails with a width of 513 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	513 mm (A-B, B-C)	1	04304	49,99



Waldner SCALA fume hood with rear wall installation - Width 1500 mm

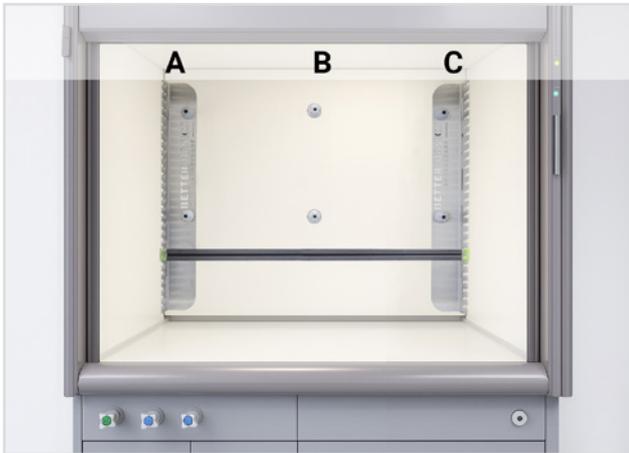


Fume hoods with a width of 1500 mm have three rows (A, B and C) of stand rod holders.

To connect the outer rows of stand rod holders (A-C), select the required number of slot rails with a width of 1401 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	1401 mm (A-C)	1	04327	79,99

Waldner SCALA fume hood with rear wall installation - Width 1500 mm

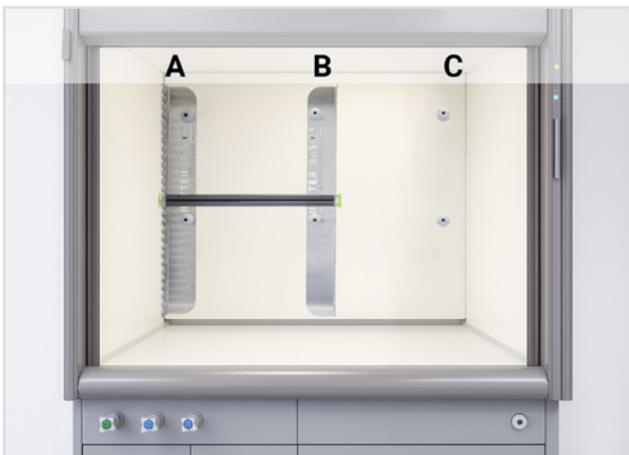

Fume hoods with a width of 1500 mm have three rows (A, B and C) of stand rod holders.


To connect the middle row of stand rod holders with one of the outer rows (A-B, B-C), select the required number of slot rails with a width of 763 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	763 mm (A-B, B-C)	1	04314	59,99



Waldner SCALA fume hood with side installation - Width 1500 mm

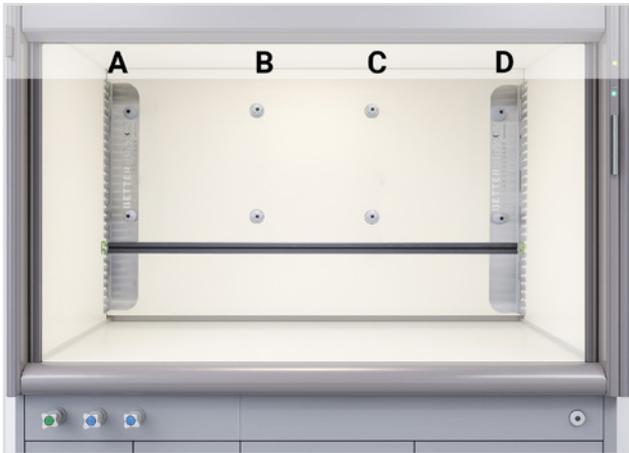

Fume hoods with a width of 1500 mm have four rows (A, B, C and D) of stand rod holders.

To connect the outer rows of stand rod holders (A-D), select the required number of slot rails with a width of 1200 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	1200 mm (A-D)	1	04322	74,99

Waldner SCALA fume hood with side installation - Width 1500 mm



Fume hoods with a width of 1500 mm have four rows (A, B, C and D) of stand rod holders.


To connect the stand rod holders in such a way that one outer row of stand rod holders remains free (A-C, B-D), select the required number of slot rails with a width of 850 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	850 mm (A-C, B-D)	1	04317	59,99



Waldner SCALA fume hood with side installation - Width 1500 mm

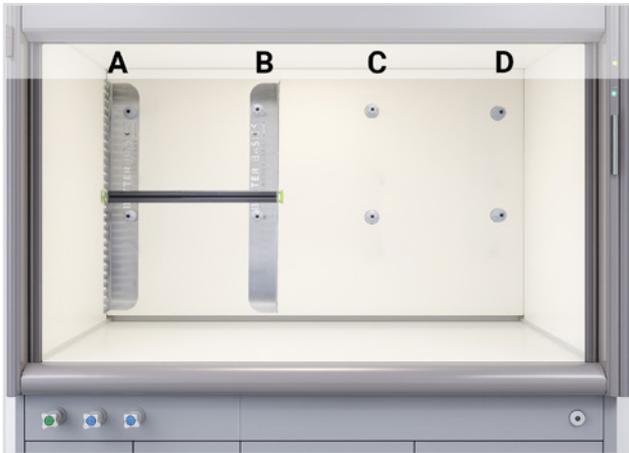


Fume hoods with a width of 1500 mm have four rows (A, B, C and D) of stand rod holders.

To connect the two left or the two right rows of stand rod holders (A-B, C-D), select the required number of slot rails with a width of 475 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	475 mm (A-B, C-D)	1	04302	49,99

Waldner SCALA fume hood with side installation - Width 1500 mm

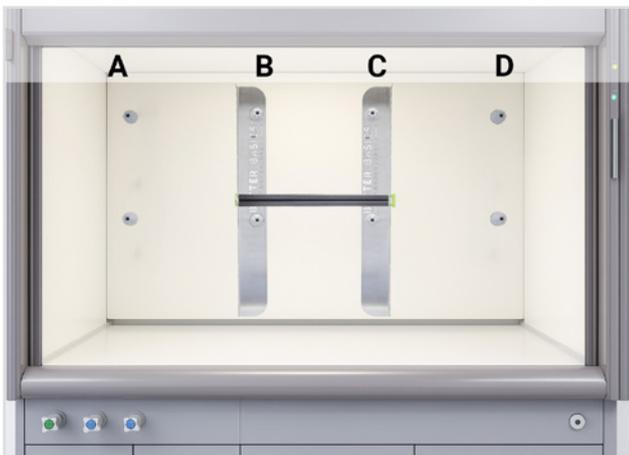

Fume hoods with a width of 1500 mm have four rows (A, B, C and D) of stand rod holders.


To connect the two middle rows of stand rod holders (B-C), select the required number of slot rails with a width of 500 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	500 mm (B-C)	1	04303	49,99



Waldner SCALA fume hood with rear wall installation - Width 1800 mm

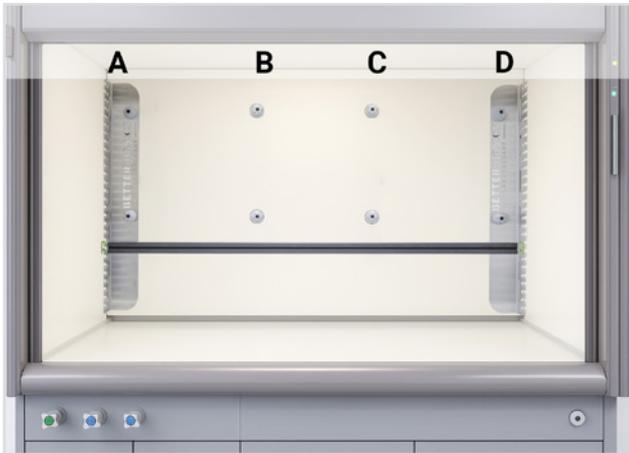

Fume hoods with a width of 1800 mm have four rows (A, B, C and D) of stand rod holders.

To connect the outer rows of stand rod holders (A-D), select the required number of slot rails with a width of 1701 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	1701 mm (A-D)	1	04330	89,99

Waldner SCALA fume hood with rear wall installation - Width 1800 mm



Fume hoods with a width of 1800 mm have four rows (A, B, C and D) of stand rod holders.


To connect the stand rod holders in such a way that one outer row of stand rod holders remains free (A-C, B-D), select the required number of slot rails with a width of 1207 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	1207 mm (A-C, B-D)	1	04324	74,99



Waldner SCALA fume hood with rear wall installation - Width 1800 mm

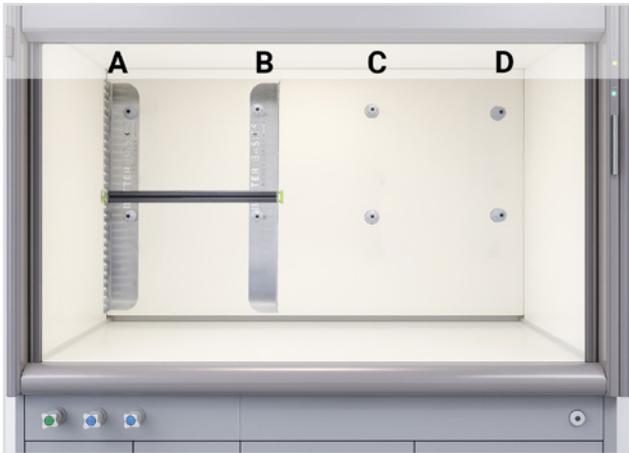


Fume hoods with a width of 1800 mm have four rows (A, B, C and D) of stand rod holders.

To connect the two left or the two right rows of stand rod holders (A-B, C-D), select the required number of slot rails with a width of 619 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	619 mm (A-B, C-D)	1	04308	54,99

Waldner SCALA fume hood with rear wall installation - Width 1800 mm

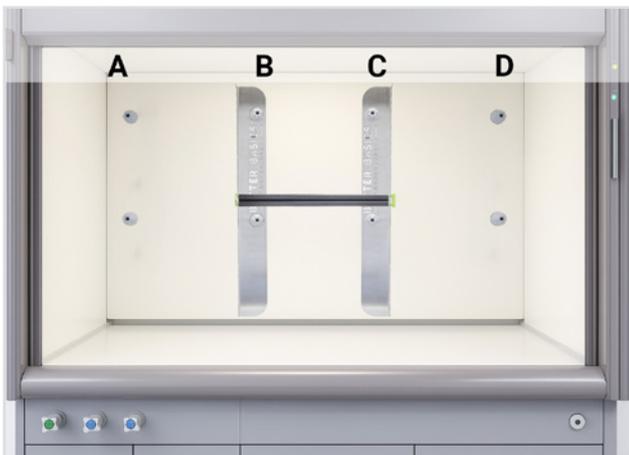

Fume hoods with a width of 1800 mm have four rows (A, B, C and D) of stand rod holders.


To connect the two middle rows of stand rod holders (B-C), select the required number of slot rails with a width of 713 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	713 mm (B-C)	1	04312	54,99



Waldner SCALA fume hood with side installation - Width 1800 mm

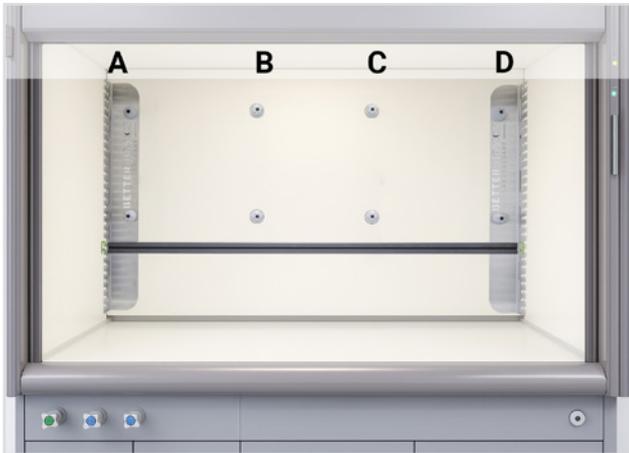

Fume hoods with a width of 1800 mm have four rows (A, B, C and D) of stand rod holders.

To connect the outer rows of stand rod holders (A-D), select the required number of slot rails with a width of 1500 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	1500 mm (A-D)	1	04328	84,99

Waldner SCALA fume hood with side installation - Width 1800 mm



Fume hoods with a width of 1800 mm have four rows (A, B, C and D) of stand rod holders.


To connect the stand rod holders in such a way that one outer row of stand rod holders remains free (A-C, B-D), select the required number of slot rails with a width of 1100 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	1100 mm (A-C, B-D)	1	04319	69,99



Waldner SCALA fume hood with side installation - Width 1800 mm

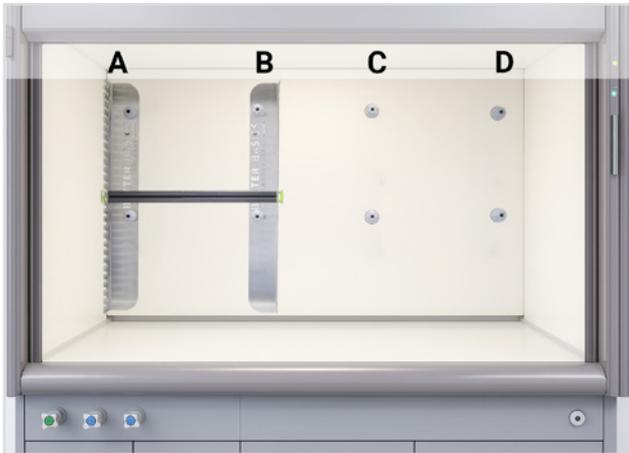


Fume hoods with a width of 1800 mm have four rows (A, B, C and D) of stand rod holders.

To connect the two left or the two right rows of stand rod holders (A-B, C-D), select the required number of slot rails with a width of 525 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	525 mm (A-B, C-D)	1	04305	49,99

Waldner SCALA fume hood with side installation - Width 1800 mm

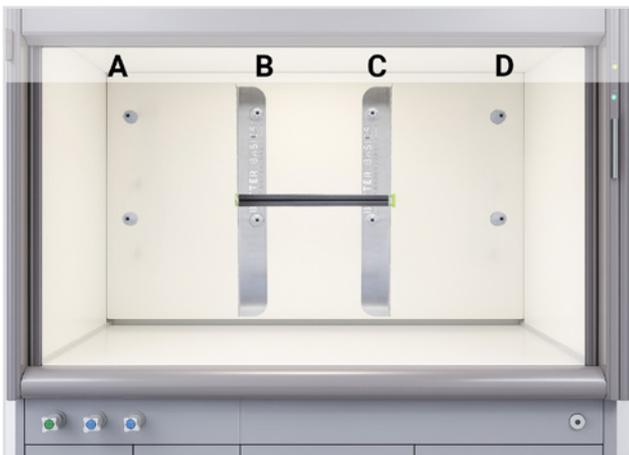

Fume hoods with a width of 1800 mm have four rows (A, B, C and D) of stand rod holders.


To connect the two middle rows of stand rod holders (B-C), select the required number of slot rails with a width of 700 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	700 mm (B-C)	1	04311	54,99



Waldner SCALA fume hood with rear wall installation - Width 2100 mm

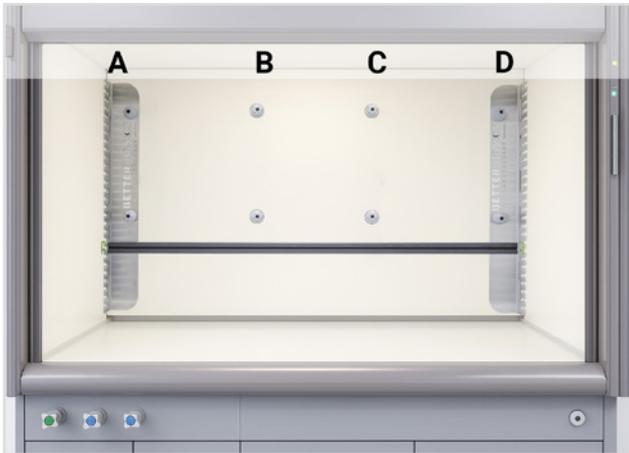

Fume hoods with a width of 2100 mm have four rows (A, B, C and D) of stand rod holders.

To connect the outer rows of stand rod holders (A-D), select the required number of slot rails with a width of 2001 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	2001 mm (A-D)	1	04332	99,99

Waldner SCALA fume hood with rear wall installation - Width 2100 mm



Fume hoods with a width of 2100 mm have four rows (A, B, C and D) of stand rod holders.


To connect the stand rod holders in such a way that one outer row of stand rod holders remains free (A-C, B-D), select the required number of slot rails with a width of 1357 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	1357 mm (A-C, B-D)	1	04326	79,99



Waldner SCALA fume hood with rear wall installation - Width 2100 mm

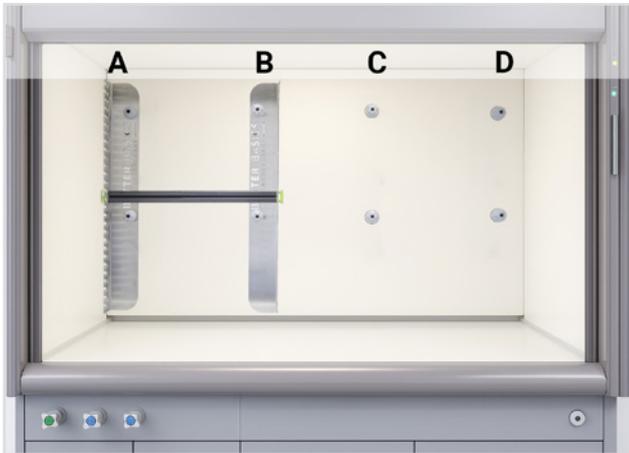


Fume hoods with a width of 2100 mm have four rows (A, B, C and D) of stand rod holders.

To connect the two left or the two right rows of stand rod holders (A-B, C-D), select the required number of slot rails with a width of 769 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	769 mm (A-B, C-D)	1	04315	59,99

Waldner SCALA fume hood with rear wall installation - Width 2100 mm

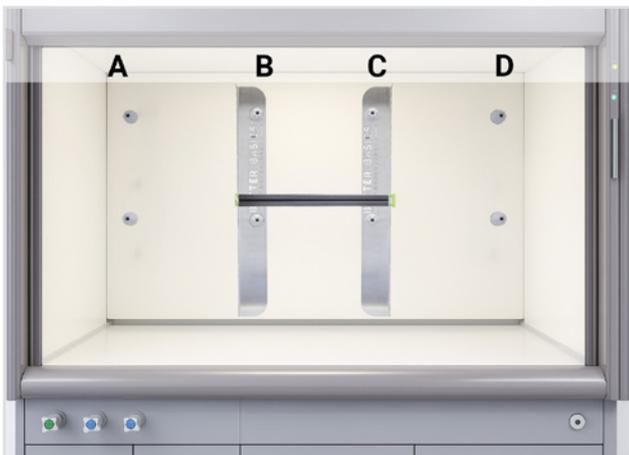

Fume hoods with a width of 2100 mm have four rows (A, B, C and D) of stand rod holders.


To connect the two middle rows of stand rod holders (B-C), select the required number of slot rails with a width of 713 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	713 mm (B-C)	1	04313	54,99



Waldner SCALA fume hood with side installation - Width 2100 mm

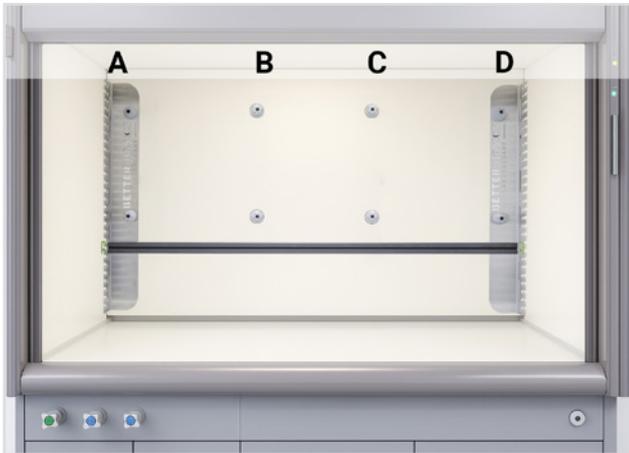

Fume hoods with a width of 2100 mm have four rows (A, B, C and D) of stand rod holders.

To connect the outer rows of stand rod holders (A-D), select the required number of slot rails with a width of 1800 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	1800 mm (A-D)	1	04331	94,99

Waldner SCALA fume hood with side installation - Width 2100 mm



Fume hoods with a width of 2100 mm have four rows (A, B, C and D) of stand rod holders.


To connect the stand rod holders in such a way that one outer row of stand rod holders remains free (A-C, B-D), select the required number of slot rails with a width of 1300 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	1300 mm (A-C, B-D)	1	04325	74,99



Waldner SCALA fume hood with side installation - Width 2100 mm

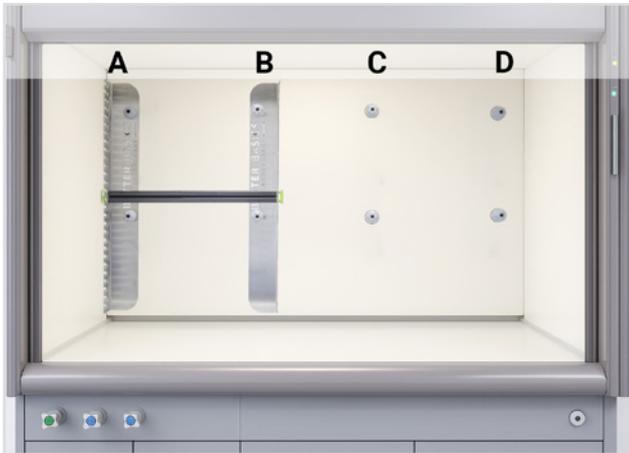


Fume hoods with a width of 2100 mm have four rows (A, B, C and D) of stand rod holders.

To connect the two left or the two right rows of stand rod holders (A-B, C-D), select the required number of slot rails with a width of 625 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	625 mm (A-B, C-D)	1	04309	54,99

Waldner SCALA fume hood with side installation - Width 2100 mm

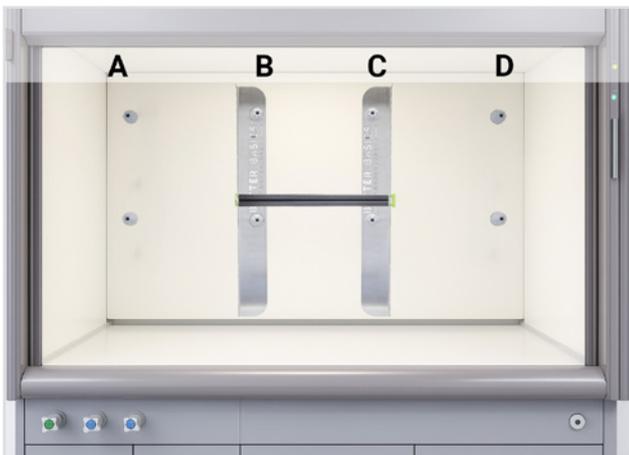

Fume hoods with a width of 2100 mm have four rows (A, B, C and D) of stand rod holders.


To connect the two middle rows of stand rod holders (B-C), select the required number of slot rails with a width of 800 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	800 mm (B-C)	1	04316	59,99



Waldner SCALA fume hood with side installation - Width 2400 mm

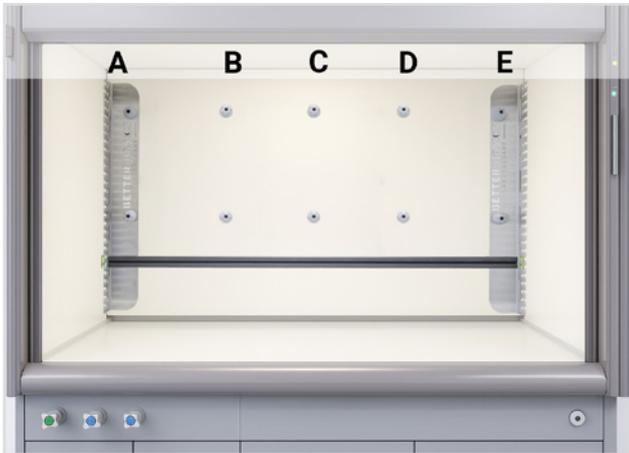


Fume hoods with a width of 2400 mm have five rows (A, B, C, D and E) of stand rod holders.

To connect the outer rows of stand rod holders (A-E), select the required number of slot rails with a width of 2101 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	2101 mm (A-E)	1	04333	99,99

Waldner SCALA fume hood with side installation - Width 2400 mm

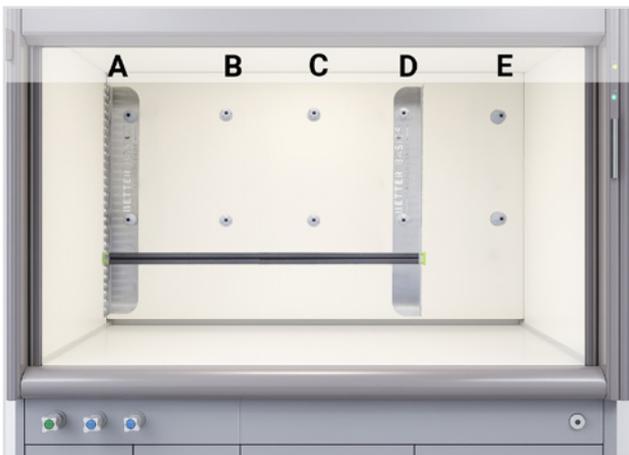

Fume hoods with a width of 2400 mm have five rows (A, B, C, D and E) of stand rod holders.


To connect the stand rod holders in such a way that one outer row of stand rod holders remains free (A-D, B-E), select the required number of slot rails with a width of 1651 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	1651 mm (A-D, B-E)	1	04329	89,99



Waldner SCALA fume hood with side installation - Width 2400 mm

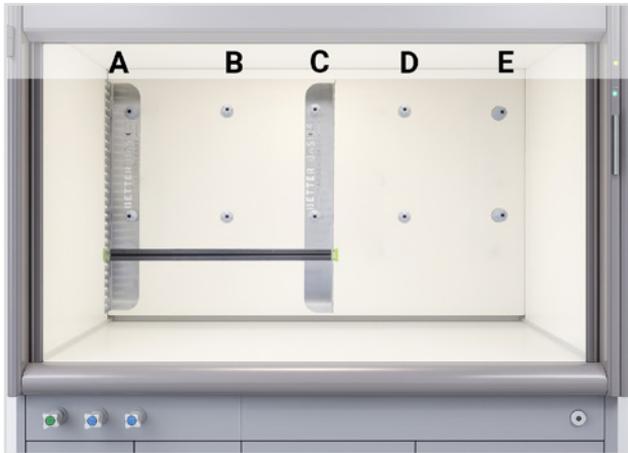


Fume hoods with a width of 2400 mm have five rows (A, B, C, D and E) of stand rod holders.

To connect the three left or three right rows of stand rod holders (A-C, C-E), select the required number of slot rails with a width of 1113 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	1113 mm (A-C, C-E)	1	04321	69,99

Waldner SCALA fume hood with side installation - Width 2400 mm

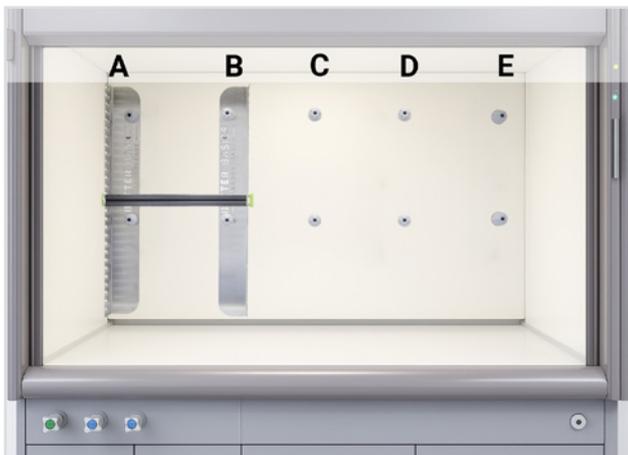

Fume hoods with a width of 2400 mm have five rows (A, B, C, D and E) of stand rod holders.


To connect the two outer left or the two outer right rows of stand rod holders (A-B, D-E), select the required number of slot rails with a width of 575 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	575 mm (A-B, D-E)	1	04306	49,99



Waldner SCALA fume hood with side installation - Width 2400 mm

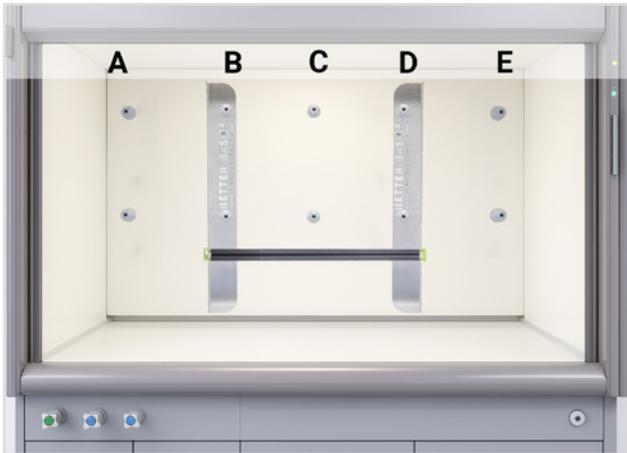


Fume hoods with a width of 2400 mm have five rows (A, B, C, D and E) of stand rod holders.

To connect the three middle rows of stand rod holders (B-D), select the required number of slot rails with a width of 1201 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	1201 mm (B-D)	1	04323	74,99

Waldner SCALA fume hood with side installation - Width 2400 mm

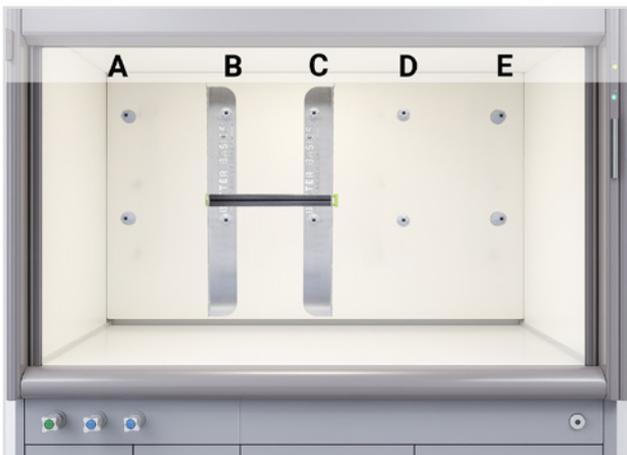

Fume hoods with a width of 2400 mm have five rows (A, B, C, D and E) of stand rod holders.

To connect the second and third or the third and fourth rows of stand rod holders (B-C, C-D), select the required number of slot rails with a width of 663 mm.

** Image for illustrative purposes only*

Item	Type	PU	Item no.	Price in EUR
SmartEdge elite WA 001	elite	1	04009	229,99
Slot rail for SmartEdge WA 001	663 mm (B-C, C-D)	1	04310	54,99



Slot rails for SmartEdge elite WA 001

The slot rail can be inserted at various positions within the SmartEdge elite WA 001. A larger number of slot rails provides more options for module placement through optimal vertical arrangement, offering greater flexibility for personalising the SmartEdge and adapting the module configuration to the individual needs of each user.

Item	Type	PU	Item no.	Price in EUR
Slot rail for SmartEdge WA 001	in RWI 1200 A-C	1	04320	69,99
Slot rail for SmartEdge WA 001	in RWI 1200 A-B or B-C	1	04307	54,99
Slot rail for SmartEdge WA 001	in SI 1200 A-C	1	04318	64,99
Slot rail for SmartEdge WA 001	in SI 1200 A-B or B-C	1	04304	49,99
Slot rail for SmartEdge WA 001	in RWI 1500 A-C	1	04327	79,99
Slot rail for SmartEdge WA 001	in RWI 1500 A-B or B-C	1	04314	59,99
Slot rail for SmartEdge WA 001	in SI 1500 A-D	1	04322	74,99
Slot rail for SmartEdge WA 001	in SI 1500 B-C	1	04303	49,99
Slot rail for SmartEdge WA 001	in SI 1500 A-B or C-D	1	04302	49,99
Slot rail for SmartEdge WA 001	in SI 1500 A-C or B-D	1	04317	59,99
Slot rail for SmartEdge WA 001	in RWI 800 A-D	1	04330	89,99
Slot rail for SmartEdge WA 001	in RWI 1800 B-C	1	04312	54,99
Slot rail for SmartEdge WA 001	in RWI 1800 A-B or C-D	1	04308	54,99
Slot rail for SmartEdge WA 001	in RWI 1800 A-C or B-D	1	04324	74,99
Slot rail for SmartEdge WA 001	in SI 1800 A-D	1	04328	84,99
Slot rail for SmartEdge WA 001	in SI 1800 B-C	1	04311	54,99
Slot rail for SmartEdge WA 001	in SI 1800 A-B or C-D	1	04305	49,99
Slot rail for SmartEdge WA 001	in SI 1800 A-C or B-D	1	04319	69,99
Slot rail for SmartEdge WA 001	in RWI 2100 A-D	1	04332	99,99
Slot rail for SmartEdge WA 001	in RWI 2100 B-C	1	04313	54,99
Slot rail for SmartEdge WA 001	in RWI 2100 A-B or C-D	1	04315	59,99
Slot rail for SmartEdge WA 001	in RWI 2100 A-C or B-D	1	04326	79,99
Slot rail for SmartEdge WA 001	in SI 2100 A-D	1	04331	94,99
Slot rail for SmartEdge WA 001	in SI 2100 B-C	1	04316	59,99
Slot rail for SmartEdge WA 001	in SI 2100 A-B or C-D	1	04309	54,99
Slot rail for SmartEdge WA 001	in SI 2100 A-C or B-D	1	04325	74,99
Slot rail for SmartEdge WA 001	in SI 2400 A-E	1	04333	99,99
Slot rail for SmartEdge WA 001	in SI 2400 B-D	1	04323	74,99
Slot rail for SmartEdge WA 001	in SI 2400 A-B or D-E	1	04306	49,99
Slot rail for SmartEdge WA 001	in SI 2400 B-C or C-D	1	04310	54,99
Slot rail for SmartEdge WA 001	in SI 2400 A-C or C-E	1	04321	69,99
Slot rail for SmartEdge WA 001	in SI 2400 A-D or B-E	1	04329	89,99



SmartFrame elite WA 001



SmartFrame elite WA 001 mounted on a lab stand system, Demonstrators, lab stand system and modules are not included in the scope of delivery.

If your laboratory fume hood is a customised design or an adaptation of the specified standard type, or if you are using another WALDNER fume hood model, please contact us. We will gladly check compatibility and develop a tailored solution for you if needed. Upon request, we can manufacture the SmartFrame to precisely fit your requirements.

The SmartFrame elite WA 001 was developed for use in the WALDNER Vision 1500 fume hood. It consists of a polished stainless steel frame. Using separately available stand locks, the SmartFrame can be moved away from the rear wall of the fume hood towards the sash on stand rods that are located in the standard stand rod holders in Vision series fume cupboards. (NOTE: Stand rods are not included in the scope of delivery).


Every third slot rail holder on both sides of the SmartFrame elite has a recess at the front to make it easier to align the slot rails horizontally. All modules are fully compatible with the SmartFrame

The SmartFrame elite WA 001 comprises the following components:

- 1 x frame
- 5 x slot rails without handles



SmartFrame elite WA 001 mounted on a lab stand system, Demonstrators, lab stand system and modules are not included in the scope of delivery.

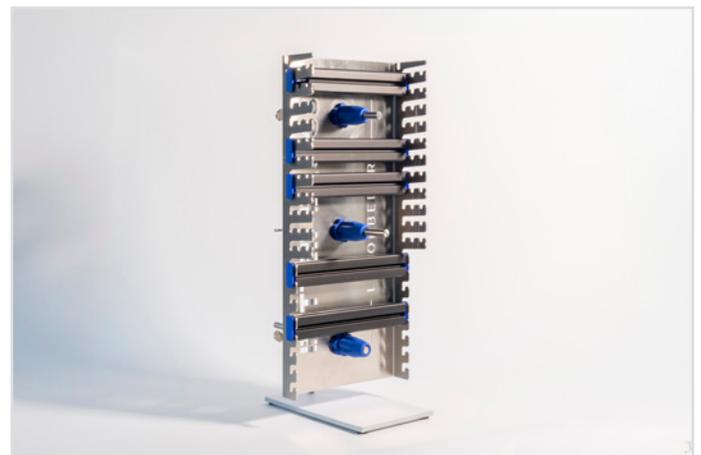

SmartFrame elite WA 001 mounted on a lab stand system, Lab stand system is not included in the scope of delivery.

Item	Type	PU	Item no.	Price in EUR
SmartFrame WA 001	elite	1	05002	579,00
Slot rail for SmartFrame WA 001		1	05302	54,99



SmartFrame elite WA 003



SmartFrame elite WA 003,
Demonstrators and modules are not included in the scope of delivery.

If your laboratory fume hood is a customised design or an adaptation of the specified standard type, or if you are using another WALDNER fume hood model, please contact us. We will gladly check compatibility and develop a tailored solution for you if needed. Upon request, we can manufacture the SmartFrame to precisely fit your requirements.



SmartFrame elite WA 003, corner view.

The SmartFrame elite WA 003 was developed for use in various WALDNER fume cupboards in the Scala series and can be installed in the following models:

- Bench-mounted fume hood / Secuflow bench-mounted fume hood 1500
- Bench-mounted fume hood / Secuflow bench-mounted fume hood 2100
- Low ceiling bench-mounted fume hood / Secuflow low ceiling bench-mounted fume hood 1500
- Low ceiling bench-mounted fume hood / Secuflow low ceiling bench-mounted fume hood 2100
- Secuflow bench-mounted fume hood with side installation 1200
- Secuflow bench-mounted fume hood with side installation 1500
- Secuflow bench-mounted fume hood with side installation 1800
- Secuflow bench-mounted fume hood with side installation 2100
- Secuflow bench-mounted fume hood with side installation 2400
- Walk-in fume hood with side installation 1200
- Walk-in fume hood with side installation 1500
- Walk-in fume hood with side installation 1800
- Walk-in fume hood with side installation 2100
- Walk-in fume hood with side installation 2400

The SmartFrame elite WA 003 consists of a polished stainless steel frame. Using separately available stand locks, the SmartFrame can be moved away from the rear wall of the fume hood towards the sash on stand rods that are located in the standard lab stand holders in fume hoods of the Scala and Vision series. (NOTE: Stand rods are not included in the scope of delivery).

Every third slot rail holder on both sides of the SmartFrame elite has a recess at the front to make it easier to align the slot rails horizontally. All modules are fully compatible with this SmartFrame.

This SmartFrame elite WA 003 comprises the following components:

- 1 x frame
- 5 x slot rails without handles

Item	Type	PU	Item no.	Price in EUR
SmartFrame WA 003	elite	1	05012	599,00
Slot rail for SmartFrame WA 003		1	05308	54,99



SmartRail WM 001



*SmartRail WM 001 mounted,
Demonstrators and modules are not included in the scope
of delivery.*

The SmartRail WM 001 is an innovative solution for the effective organization of laboratory workspaces. This specially designed slot rail is made of black anodised aluminium and is intended for secure mounting on walls or vertical surfaces using screws. At the ends (or possibly in the middle of the rail), there are screw mounts made of PA12 and possibly stainless steel, through which two screws are used to ensure stable fixation of the rail.

SmartRail WM 001 is supplied with fixing materials (screws and dowels), which enable it to be mounted on walls made of masonry, concrete, wood and similar materials.

Components of the SmartRail WM 001:

- Anodised aluminium slot rail
- Two or more sets of mounting materials (consisting of PA12 and possibly stainless steel screw mounts, as well as appropriate mounting hardware) for installation

The SmartRail WM 001 is compatible with all modules, allowing for flexible and customizable placement of these within the slot rail. This solution not only frees up valuable space in the laboratory but also significantly enhances organization and efficiency.

Item	Type	PU	Item no.	Price in EUR
SmartRail WM 001	300 mm	1	03048	79,99
SmartRail WM 001	600 mm	1	03049	99,99
SmartRail WM 001	900 mm	1	03050	119,99
SmartRail WM 001	1200 mm	1	03051	139,99
SmartRail WM 001	1500 mm	1	03052	159,99



SmartRail MB 001



SmartRail MB 001 in a fume hood, mounted on a lab stand system, Demonstrators, lab stand system and modules are not included in the scope of delivery.

The SmartRail MB 001 has been specifically developed for use with laboratory furniture featuring stand rod systems of 12 mm or 13 mm. The design consists of a black anodised aluminium slot rail and side brackets made of PA12 and stainless steel, equipped with a sturdy screw clamp to ensure stable mounting on the support system. The versatile SmartRail offers flexible placement options in the laboratory. A larger number of SmartRail provides various module placement options and allows the laboratory environment to be tailored to specific requirements. All modules are fully compatible with SmartRail.

The SmartRail MB 001 can be used on lab stand systems both in the fume hood and on the workbench.



SmartRail MB 001 mounted, without equipment, lab stand system is not included in the scope of delivery.



SmartRail MB 001, corner view.

Item	Type	PU	Item no.	Price in EUR
SmartRail MB 001	300 mm	1	03023	79,99
SmartRail MB 001	600 mm	1	03024	99,99
SmartRail MB 001	900 mm	1	03025	119,99
SmartRail MB 001	1200 mm	1	03026	139,99
SmartRail MB 001	1500 mm	1	03027	159,99



SmartFrame elite MB 001




SmartFrame elite MB 001, front view, Demonstrators, lab stand system and modules are not included in the scope of delivery.

The SmartFrame elite MB 001 is designed for use in fume hoods. It consists of a polished stainless steel frame. The SmartFrame can be mounted on standard 12mm or 13mm lab stand system using PA12 brackets on the rear. Depending on the application, it can be mounted close to the back of the fume hood, free-standing in the fume hood or on a laboratory bench using a lab stand system. (NOTE: lab stand rods are not supplied) Every third slot rail holder on both sides of the SmartFrame elite has a recess at the front to make it easier to align the slot rails horizontally. All modules are fully compatible with this SmartFrame.


The SmartFrame elite MB 001 consists of the following components:

- 1 x frame
- 3 x slot rails with handles
- 3 x PA12 universal support rod brackets

The SmartFrame elite MB 001 can be used on lab stand systems both in the fume hood and on the workbench.

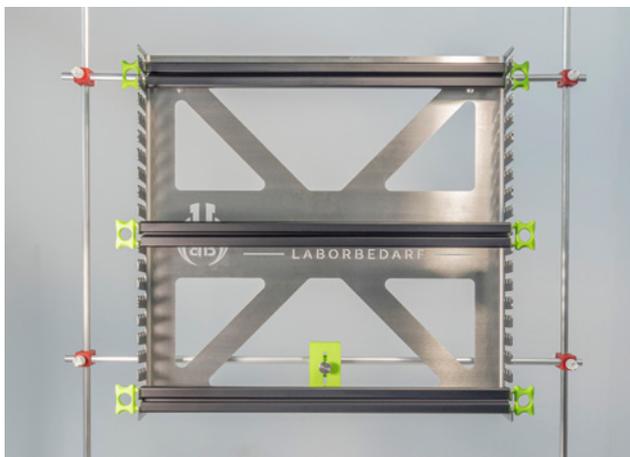

SmartFrame elite MB 001 mounted on a lab stand system, Lab stand system is not included in the scope of delivery.

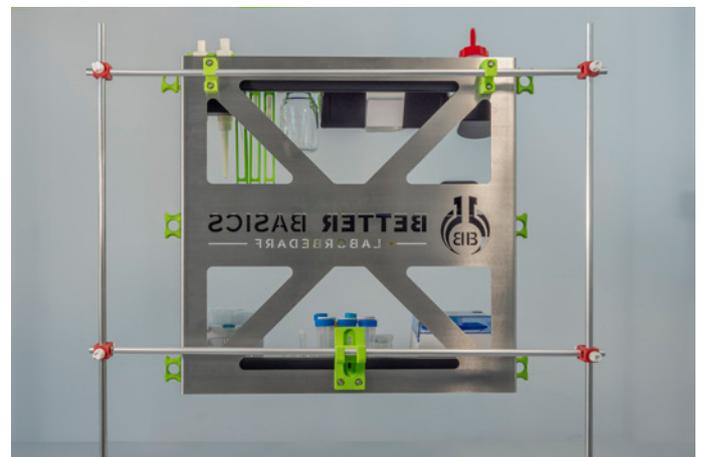

SmartFrame elite MB 001, view from behind with demonstrators and modules, Demonstrators, lab stand system and modules are not included in the scope of delivery.

Item	Type	PU	Item no.	Price in EUR
SmartFrame elite MB 001	635 x 615 x 120 mm (WxHxD)	1	05005	599,00
SmartFrame elite mini MB 001	450 x 380 x 120 mm (WxHxD)	1	05006	469,00
Slot rail for SmartFrame MB 001		1	05312	74,99
Slot rail for SmartFrame mini MB 001		1	05313	54,99



SmartFrame elite VESA 100



SmartFrame elite VESA 100, front view.

The SmartFrame elite Vesa 100 allows you to mount a SmartFrame elite on a VESA 100 mount. It consists of a polished stainless steel frame. The SmartFrame elite can be attached to the VESA 100 mount using screws. (NOTE: VESA 100 mount not included) Every third slot rail bracket on either side of the SmartFrame elite has a recess at the front to make it easier to align the slot rails horizontally. All modules are fully compatible with SmartFrame.

The SmartFrame elite VESA 100 consists of the following components:

- 1 x frame
- 3 x slot rails with handles

The SmartFrame elite VESA 100 can be used on VESA 100 mounts both in the fume hood and on the workbench.

Item	Type	PU	Item no.	Price in EUR
SmartFrame elite VESA 100	635 x 615 x 120 mm (WxHxD)	1	05001	549,00
Slot rail for SmartFrame VESA 100		1	05311	74,99

Additional Accessories

SmartLock



The SmartLock is used to additionally secure a SmartRail WA 001 to prevent unintentional disassembly and to ensure stable mounting on the WALDNER functional rail.

Item	Type	PU	Item no.	Price in EUR
SmartLock		1	09004	4,99

Support rod lock for SmartFrame WA 001 and WA 003



The support rod lock for SmartFrame WA 001 and WA 003 allows the SmartFrames to be mounted on support rods that are inserted into the support rod holders of the corresponding fume cupboards. The SmartFrames can therefore be positioned closer to the user.

Item	Type	PU	Item no.	Price in EUR
Support rod lock for SmartFrame WA 001 and WA 003		1	05307	16,99

Customising

Better Basics Laborbedarf offers an unparalleled approach to laboratory workstation organisation with its customisable and modular solutions. Below is an overview of how our customisation services span both the SmartRack® and SmartIntegrate product lines:

1. Customisation of the SmartRack®

For the SmartRack®

Our SmartRack® system is fully customisable to meet the specific requirements of your workspace. The dimensions of the SmartRack® can be individually adapted to the size of your equipment and workflows, ensuring seamless integration into your laboratory workstation. At the same time, the modular compatibility and reliability of your tailor-made SmartRack® are guaranteed through appropriately dimensioned slot rails. From compact configurations to larger frames, we enhance the organisation and efficiency of your laboratory.

2. Customisation of SmartIntegrate Products

a) Bespoke SmartFrame for Fume Hoods

The SmartFrame optimises the organisation within fume hoods, offering a structured and safe arrangement for laboratory glassware, tools, and equipment. Custom-built SmartFrames provide a perfect fit for your specific workspace.

b) Tailored SmartRail for Lab Benches

SmartRails bring adaptability and precision to your laboratory benches. Designed for both high-throughput workflows and compact setups, they ensure reduced clutter and enhanced efficiency. Each SmartRail is tailored to your preferred length and mounting specifications.

c) Custom Solutions for Safety Cabinets

For safety cabinets, our custom solutions include integrated SmartRails, SmartFrames and SmartEdges designed to maximise space utilisation while supporting your workflows. Each component is built to ensure easy cleaning and adherence to safety regulations.

3. Custom Modules

Tailored for any Laboratory Equipment

Using advanced 3D printing technology, we create modules and holders in any shape or size to fit your laboratory equipment. From pipette holders to mounts for specialised glassware, our bespoke designs cater to your unique requirements. Our rapid prototyping process ensures swift delivery of tailored solutions.

Why Choose Better Basics for Customisation?

- 1. Precision Craftsmanship:** Each component is tailored to your exact specifications, guaranteeing optimal performance and seamless integration.
- 2. Fast Delivery:** Our in-house design and manufacturing teams ensure swift turnaround from concept to delivery.
- 3. Transparent Pricing:** Costs for customisation are agreed upon in advance, providing clarity and control.
- 4. Scalability:** Whether expanding an existing laboratory or building a new one, our products adapt to meet your growing needs.

Transform Your Laboratory Workspace

With Better Basics Laborbedarf, your laboratory can benefit from customised SmartRack® or SmartIntegrate products.

Visit customising.better-basics-laborbedarf.de or contact our team today to unlock the full potential of your laboratory space.

This concise yet comprehensive approach ensures your laboratory benefits from enhanced efficiency, organisation, and safety—all tailored to your specific needs.

Examples for realised customisations

Special frame sizes

SmartRack® frame - customisable in width as well as height



SmartRack® frame - extra wide

For a special requirement we manufactured an extra wide frame of the SmartRack®. Due to the width, this SmartRack® makes optimal use of the working space in our customer's laboratory.

Special module sizes

Holder for multiple glassware



Double holder for glassware with ST10

For a project, a customer needed a double holder for glassware ST10. The module allows glassware with standard ground joint ST10 to be easily hung next to each other. The specially shaped recess prevents the glassware from falling out.

Of course, adaptations of this kind are also conceivable for other modules, e.g. storage boxes in individual sizes.

Special laboratory equipment

Dispenser for cylindrical Petri dish multi-packs



For a globally operating biotechnology company, a dispenser for cylindrical Petri dish multi-packs was developed.

The dispenser not only enables easy removal of the packs but also provides a quick overview of the remaining stock.

Item Overview

Waldner SmartRail (work bench)

Item	Type	Item no.	PU	Price in EUR	Page
SmartRail WA 001	300 mm	03002	1	69,99	7
SmartRail WA 001	600 mm	03003	1	89,99	7
SmartRail WA 001	900 mm	03004	1	109,99	7
SmartRail WA 001	1200 mm	03005	1	129,99	7
SmartRail WA 001	1500 mm	03006	1	149,99	7
SmartRail WA 002	300 mm	03007	1	79,99	8
SmartRail WA 002	600 mm	03008	1	99,99	8
SmartRail WA 002	900 mm	03009	1	119,99	8
SmartRail WA 002	1200 mm	03010	1	139,99	8
SmartRail WA 002	1500 mm	03011	1	159,99	8
SmartRail WA 003	300 mm	03012	1	79,99	9
SmartRail WA 003	600 mm	03013	1	99,99	9
SmartRail WA 003	900 mm	03014	1	119,99	9
SmartRail WA 003	1200 mm	03015	1	139,99	9
SmartRail WA 003	1500 mm	03016	1	159,99	9

Waldner SmartEdge (fume hood)

Item	Type	Item no.	PU	Price in EUR	Page
Slot rail for SmartEdge WA 001	in RWI 1200 A-C	04320	1	69,99	11
Slot rail for SmartEdge WA 001	in RWI 1200 A-B or B-C	04307	1	54,99	11
Slot rail for SmartEdge WA 001	in SI 1200 A-C	04318	1	64,99	12
Slot rail for SmartEdge WA 001	in SI 1200 A-B or B-C	04304	1	49,99	12
Slot rail for SmartEdge WA 001	in RWI 1500 A-C	04327	1	79,99	13
Slot rail for SmartEdge WA 001	in RWI 1500 A-B or B-C	04314	1	59,99	13
Slot rail for SmartEdge WA 001	in SI 1500 A-D	04322	1	74,99	14
Slot rail for SmartEdge WA 001	in SI 1500 B-C	04303	1	49,99	15
Slot rail for SmartEdge WA 001	in SI 1500 A-B or C-D	04302	1	49,99	15
Slot rail for SmartEdge WA 001	in SI 1500 A-C or B-D	04317	1	59,99	14
Slot rail for SmartEdge WA 001	in RWI 800 A-D	04330	1	89,99	16
Slot rail for SmartEdge WA 001	in RWI 1800 B-C	04312	1	54,99	17
Slot rail for SmartEdge WA 001	in RWI 1800 A-B or C-D	04308	1	54,99	17
Slot rail for SmartEdge WA 001	in RWI 1800 A-C or B-D	04324	1	74,99	16
Slot rail for SmartEdge WA 001	in SI 1800 A-D	04328	1	84,99	18
Slot rail for SmartEdge WA 001	in SI 1800 B-C	04311	1	54,99	19
Slot rail for SmartEdge WA 001	in SI 1800 A-B or C-D	04305	1	49,99	19
Slot rail for SmartEdge WA 001	in SI 1800 A-C or B-D	04319	1	69,99	18
Slot rail for SmartEdge WA 001	in RWI 2100 A-D	04332	1	99,99	20
Slot rail for SmartEdge WA 001	in RWI 2100 B-C	04313	1	54,99	21
Slot rail for SmartEdge WA 001	in RWI 2100 A-B or C-D	04315	1	59,99	21

Item	Type	Item no.	PU	Price in EUR	Page
Slot rail for SmartEdge WA 001	in RWI 2100 A-C or B-D	04326	1	79,99	20
Slot rail for SmartEdge WA 001	in SI 2100 A-D	04331	1	94,99	22
Slot rail for SmartEdge WA 001	in SI 2100 B-C	04316	1	59,99	23
Slot rail for SmartEdge WA 001	in SI 2100 A-B or C-D	04309	1	54,99	23
Slot rail for SmartEdge WA 001	in SI 2100 A-C or B-D	04325	1	74,99	22
Slot rail for SmartEdge WA 001	in SI 2400 A-E	04333	1	99,99	24
Slot rail for SmartEdge WA 001	in SI 2400 B-D	04323	1	74,99	26
Slot rail for SmartEdge WA 001	in SI 2400 A-B or D-E	04306	1	49,99	25
Slot rail for SmartEdge WA 001	in SI 2400 B-C or C-D	04310	1	54,99	26
Slot rail for SmartEdge WA 001	in SI 2400 A-C or C-E	04321	1	69,99	25
Slot rail for SmartEdge WA 001	in SI 2400 A-D or B-E	04329	1	89,99	24
SmartEdge elite WA 001	elite	04009	1	229,99	10

Waldner SmartFrame (fume hood)

Item	Type	Item no.	PU	Price in EUR	Page
Slot rail for SmartFrame WA 001		05302	1	54,99	28
Slot rail for SmartFrame WA 003		05308	1	54,99	29
SmartFrame WA 001	elite	05002	1	579,00	28
SmartFrame WA 003	elite	05012	1	599,00	29

Universal Solutions (work bench/fume hood)

Item	Type	Item no.	PU	Price in EUR	Page
Slot rail for SmartFrame MB 001		05312	1	74,99	32
Slot rail for SmartFrame mini MB 001		05313	1	54,99	32
Slot rail for SmartFrame VESA 100		05311	1	74,99	33
SmartFrame elite MB 001	635 x 615 x 120 mm (WxHxD)	05005	1	599,00	32
SmartFrame elite mini MB 001	450 x 380 x 120 mm (WxHxD)	05006	1	469,00	32
SmartFrame elite VESA 100	635 x 615 x 120 mm (WxHxD)	05001	1	549,00	33
SmartRail MB 001	300 mm	03023	1	79,99	31
SmartRail MB 001	600 mm	03024	1	99,99	31
SmartRail MB 001	900 mm	03025	1	119,99	31
SmartRail MB 001	1200 mm	03026	1	139,99	31
SmartRail MB 001	1500 mm	03027	1	159,99	31
SmartRail WM 001	300 mm	03048	1	79,99	30
SmartRail WM 001	600 mm	03049	1	99,99	30
SmartRail WM 001	900 mm	03050	1	119,99	30
SmartRail WM 001	1200 mm	03051	1	139,99	30
SmartRail WM 001	1500 mm	03052	1	159,99	30

Additional Accessories

Item	Type	Item no.	PU	Price in EUR	Page
SmartLock		09004	1	4,99	34
Support rod lock for SmartFrame WA 001 and WA 003		05307	1	16,99	34

Net prices plus statutory VAT. | Valid from 01.01.2026

The general terms and conditions of Better Basics Laborbedarf GmbH apply.

Technical Information

Technical information:

- All products are washable with water and mild detergent.
- In addition all components can be wiped with alcohol-based disinfectants.
- The plastics and metals used are 100% recyclable.
- Our products are autoclavable up to 134 °C and 3 bar.
- The plastic elements are resistant to heat deformation under normal laboratory conditions.
- To gently clean the anodised aluminium rail, we recommend using a pH-neutral hand-washing detergent. Ensure that the hand wash detergent is odourless and contains no perfumes, greases or additives. Clean the slot rail by hand and rinse thoroughly with clean water.

Note on joint tolerances for ST14:

- Maximum outside diameter at the narrowest point of the cone joint: 17.2 mm (larger diameters are too wide for the holder and cannot be held)
- Minimum outer diameter at the widest point of the cone joint: 18.0 mm (smaller diameters are too narrow for the holder and would fall out)

Note on joint tolerances for ST29:

- Maximum outside diameter at the narrowest point of the cone joint: 31.4 mm (larger diameters are too wide for the holder and cannot be held)
- Minimum outer diameter at the widest point of the cone joint: 32.8 mm (smaller diameters are too narrow for the holder and would fall out)

IMPORTANT TO NOTE: The products of the SmartIntegrate line, the SmartRack® and SmartRack® mini as well as the SmartCase and SmartStand are not suitable for the permanent storage of hazardous substances.

We reserve the right to make technical changes, material or model changes, errors and misprints excepted. If you have any further questions, please do not hesitate to contact us at contact@better-basics.de.

For more detailed information, please visit our website:

Functionality of SLS 3D printing:	www.sls-3d-printing.better-basics-laborbedarf.de
Resistance to UV light:	www.uv-resistance.better-basics-laborbedarf.de
Resistance to aggressive media:	www.material-resistance.better-basics-laborbedarf.de
Sustainability and design-efficient manufacturing:	www.sustainability.better-basics-laborbedarf.de



Photo: Building of the Better Basics Laborbedarf GmbH in Dresden

YOU ARE INTERESTED IN OUR PRODUCTS OR YOU HAVE FURTHER QUESTIONS?

Get in touch with us:

contact@better-basics.de

You can also find more information and news about our product innovations at

www.Better-Basics-Laborbedarf.de

Better Basics Laborbedarf GmbH

Loebtauer Str. 69

01159 Dresden

Germany

T: +49 (0) 178 9617 577 (Direct customers)

T: +49 (0) 176 6233 8026 (Sales partners)

E: contact@better-basics.de

W: www.Better-Basics-Laborbedarf.de

Our request to you

Our products such as our SmartRack® and our SmartIntegrate products are only as good as their users. With our roots in the Leibniz Institute for Polymer Research Dresden, our company Better Basics Laborbedarf stands for the tradition of "Made in Germany" quality. We feel in multiple aspects committed to this quality label for German companies, which stands for high innovative strength. Even beyond our product catalogue, the same applies for us at Better Basics: Just contact us if you have ideas for extensions or improvements or if you notice something in your daily laboratory routine that we should optimise right away.

We want to create the best possible working environment for you in the laboratory and we are grateful for any comments and inspiration. We can jointly shape the future of laboratory work by cooperating with scientists and researchers. Just like you, we're always looking for a way to make our world a little better every day.

Impress:

Better Basics Laborbedarf GmbH
Loebtauer Str. 69, G-01159 Dresden

E-Mail: info@better-basics.de

CEO: Mario Schneider

Commercial register: HRB 39524; Register court: Amtsgericht Dresden

VAT registration number: DE327508812

