

SmartIntegrate

Integrative organisational solutions from
Better Basics for Momoline laboratory equipment

SmartEdge - Fume hood

SmartRail - Laboratory workbench

Contact us for a personalised consultation:
Contact@Better-Basics.de



SmartRail MO 002
on a Momoline laboratory bench



SmartRail MO 001

SmartRail MO 002

The SmartRail MO 001 and 002 create additional storage surfaces to free up workspace, thereby increasing efficiency and safety at the laboratory workstation.

SmartEdge

Solutions for the fume hood




SmartEdge elite MO 002 on a Momoline laboratory fumehood

The SmartEdge elite MO 002 for fume hoods of the Momoline AK series consists of two stainless steel L-profiles, which are mounted on the left and right between the support stand holders and the rear wall of the fume hood. A variable number of slot rails can be inserted into these L-profiles. This creates additional storage space and a stable, adaptable organisational base that keeps laboratory materials safe, within easy reach, and clearly arranged inside the fume hood. A larger number of slot rails provides a wide range of module placement options and enables the laboratory environment to be individually adapted to specific requirements. All Better Basics modules are fully compatible with the SmartEdge.


SmartRail

Solutions for the laboratory workbench

The SmartRail MO 001 and SmartRail MO 002 have been specifically developed for use on trays from Momoline. The SmartRail MO 001 can be mounted on storage trays, while the SmartRail MO 002 is designed for mounting on reagent trays. The SmartRail is available in different lengths. Its construction consists of a black anodised aluminium slot rail with side brackets made of PA12 and stainless steel. The side connectors of the SmartRail to the trays are each designed to ensure secure and stable installation within a short time. The versatile SmartRail provides flexible placement options in the laboratory. A larger number of SmartRails offers a wide range of module placement options and allows the laboratory environment to be individually adapted to specific requirements. All Better Basics modules are fully compatible with the SmartRail.

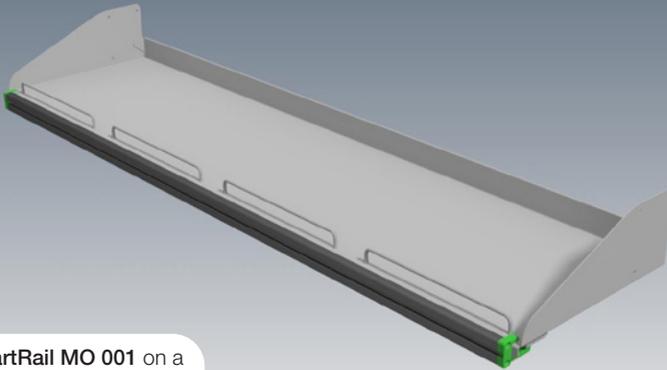

SmartRail MO 001 on a Momoline laboratory workbench

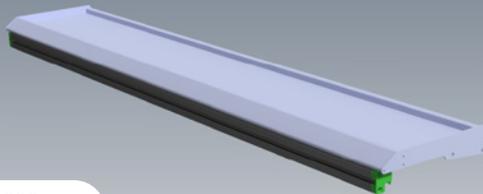

SmartRail MO 002 on a Momoline laboratory workbench

Over 100 modules

One system for endless possibilities to optimise every workplace and workflow



All Better Basics modules that can be directly mounted onto slot rails feature an integrated click mechanism, providing users with tactile feedback when positioning the holding modules within the slot rail. Since all holding modules can be freely positioned along a slot rail, you have complete flexibility in arranging them to suit your specific workflow. Once clicked into place, the modules remain securely positioned, even when the slot rail itself is moved to another working level. This ensures that both you and your work are protected against sample loss and workplace accidents.

Better Basics Laborbedarf - Unlock Your Lab Potential



Are you interested in our products or do you have further questions?

Get in touch with us: Contact@Better-Basics.de

You can also find more information and news about our product innovations at: www.Better-Basics-Laborbedarf.de

Better Basics Laborbedarf GmbH | Loebtauer Str. 69
01159 Dresden | Germany

T: +49 (0) 176 6233 8026 (Sales partner)

T: +49 (0) 178 9617 577 (direct customers)

Better Basics Laborbedarf

is a leading provider of innovative organisational systems such as the SmartRack® and the SmartIntegrate product line, revolutionising laboratory workstations. Our award-winning solutions optimise efficiency, safety, and ergonomics. By improving space utilisation and offering modular adaptability, they reduce walking distances and increase satisfaction in the laboratory. As a reliable partner, we develop customised systems tailored to individual requirements – flexible and adaptable, wherever you need them.

Impressum:

Better Basics Laborbedarf GmbH | Löbtauer Str. 69 | G-01159 Dresden
E-Mail: Info@Better-Basics.de | Geschäftsführer: Mario Schneider | Handelsregister: HRB 39524
Registergericht: Amtsgericht Dresden | Ust.-Id.-Nr.: DE327508812